




LIQUEEL®

SILICONE FOAM DRESSING WITHOUT BORDER


LIQUEEL®

Silicone foam dressing without border

LIQUEEL® consists of a waterproof, vapor-permeable outer layer made of polyurethane film [1] and a polyurethane foam layer [2] with high absorption capacity for wounds with low to heavy exudate. The soft wound-contact layer [3] helps maintain a moist wound environment to support the natural healing process. LIQUEEL® can be easily adapted to flat, convex, or concave wound sites. If necessary, foam segments [4] can be removed without the use of sterile instruments.

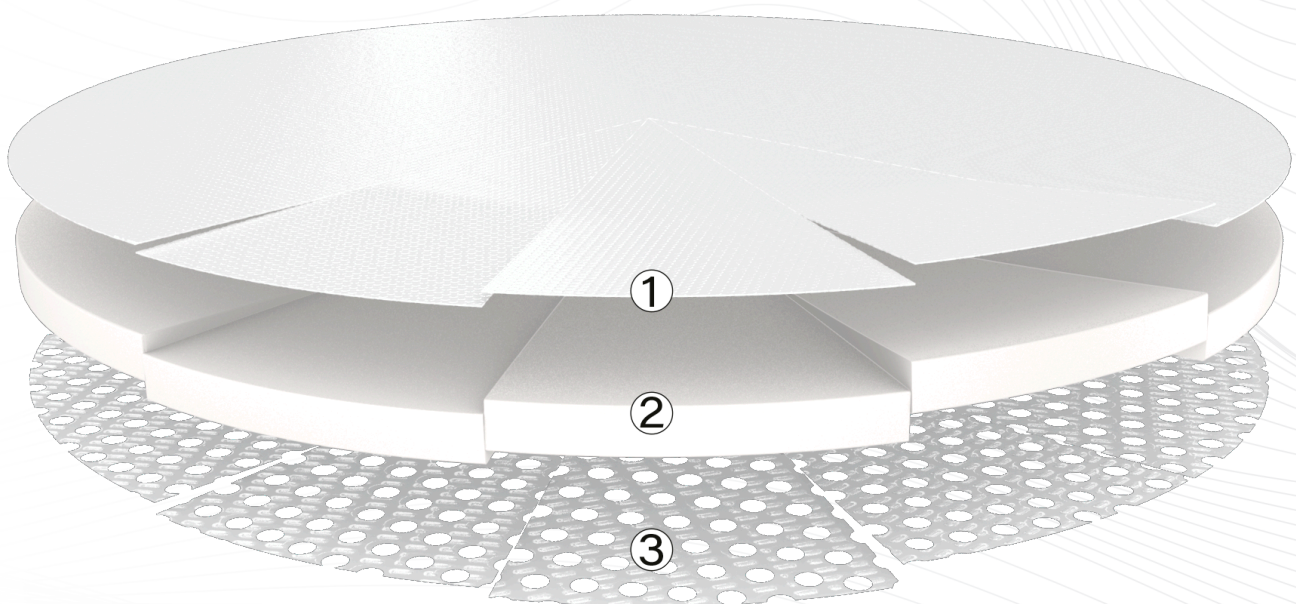**[1] LIQUEEL®**


Close-up view of the waterproof outer layer

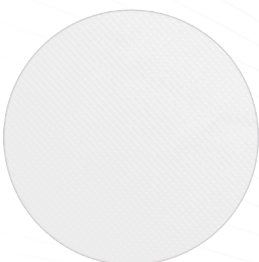**[2] LIQUEEL®**


Close-up view of the absorbent foam

**[3] LIQUEEL®**

Close-up view of the gentle-adhering silicone wound contact layer

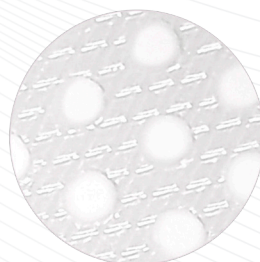**[4] LIQUEEL®**

Close-up view of the removable foam segments

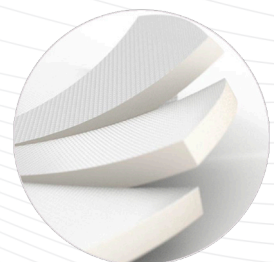

Performance

- ✓ Gentle and highly absorbent
- ✓ Minimizes wound maceration
- ✓ Atraumatic dressing changes
- ✓ Highly adaptable
- ✓ Adjustable without tools

Amount of exudate (E)

none / low 
moderate 
high 
very high 

Wound Phase (P)

Necrosis 
Fibrin coating 
Granulation 
Epithelization 
Infection 

Wound Depth (W)

shallow 
deep 
cavity 


LIQUEEL®

Smart Design for Flexible Applications

LIQUEEL® is a sterile, absorbent, self-adhesive, soft wound dressing made of polyurethane foam with a silicone coating. LIQUEEL®'s intelligent design allows the user to adapt the dressing to the specific wound condition. If necessary, individual foam segments can be removed without the use of sterile instruments.

LIQUEEL® is indicated for the treatment of wounds with mild to heavy exudate and, in the presence of exudate, helps maintain a moist wound environment to support the natural healing process.

LIQUEEL® is particularly suitable for managing convex and concave wounds, as well as wounds in specific body areas and on specific body shapes (e.g., heel, female breast, elbow, sacral region, knee); it has excellent adhesion properties and allows for atraumatic dressing changes.

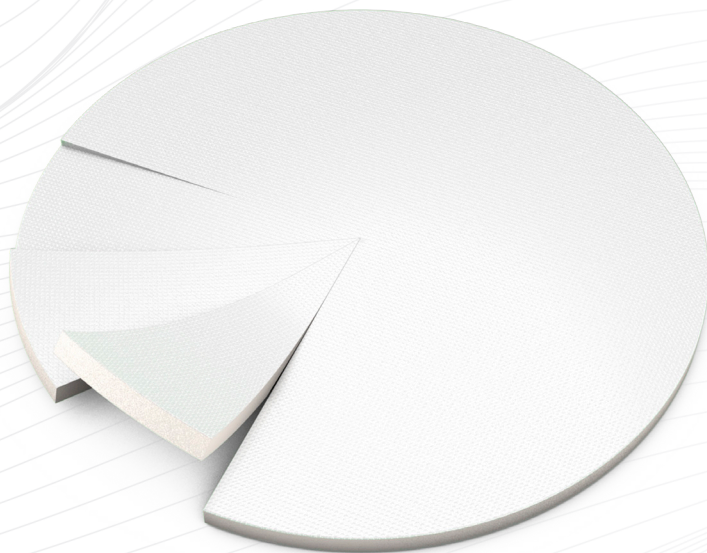

Performance

- ✓ Gentle and highly absorbent
- ✓ Minimizes wound maceration
- ✓ Atraumatic dressing changes
- ✓ Highly adaptable
- ✓ Adjustable without tools

Indications

- ✓ Pressure ulcers
- ✓ Venous leg ulcers
- ✓ Arterial leg ulcers
- ✓ Uninfected diabetic foot ulcers
- ✓ Scalds
- ✓ First- and second-degree burns
- ✓ Abrasions, lacerations, and minor cuts

	Size	Additional Information	Unit*	PZN**	REF
S	Ø 75 mm	Ø 56,5 mm bis 75 mm	1 x 10	19295624	LQ-75
M	Ø 150 mm	Ø 112,5 mm bis 150 mm	1 x 10	19295630	LQ-150
L	Ø 220 mm	Ø 165 mm bis 220 mm	1 x 10	19295647	LQ-220

LIQUEEL®

Illustration of Possible Applications


Illustration of the flexible application and optimal positioning on concave or convex areas of the body. In addition, LIQUEEL® can also be used on all other sites typically associated with chronic wounds.



Schulter



Knie

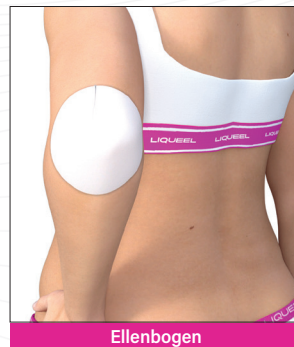

Ellenbogen



Sacral /Gluteal



Mamma



Ohr



Hand



Ferse

LIQUEEL® in Practice

Application Video


In this practical application video, we'll show you how to apply and secure LIQUEEL® to a convex body part (the knee). Scan the QR code to watch the video.



LIQUEEL®

Silicone foam dressing without border

LIQUEEL® is unique in both design and ease of use, combining many different wound dressing sizes and special shapes into a single wound dressing. LIQUEEL® really comes into its own when dealing with challenging wound locations.

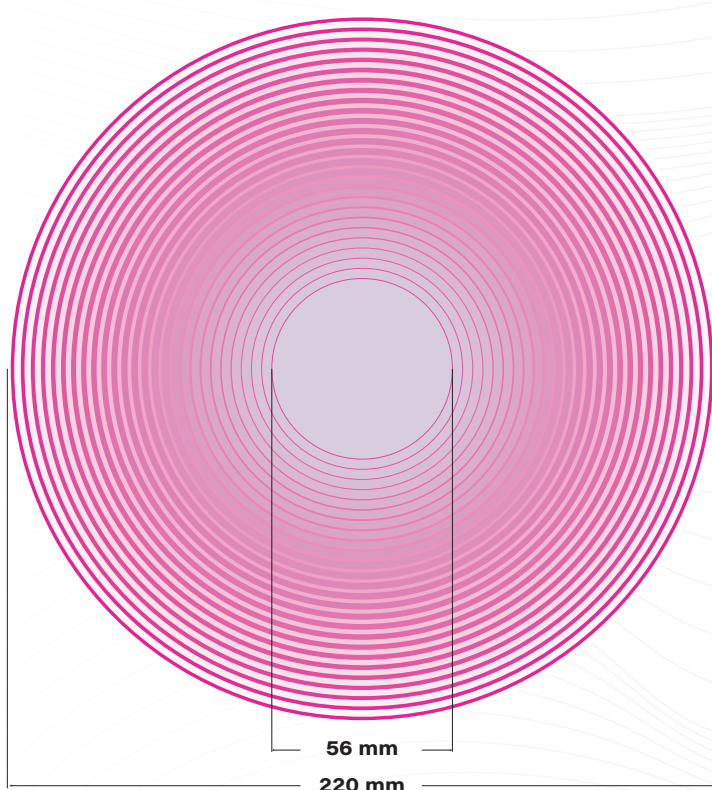

THE MULTITALENTED ONE **THREE GEOMETRIES**

LIQUEEL® is an innovative solution for the management of low- to moderately exuding wounds and can be flexibly adjusted from 56 mm to 220 mm without the use of sterile instruments.

LIQUEEL® is particularly well-suited for managing flat, convex, and concave wounds, as well as wounds on areas of the body with unique contours.

Its innovative design and wide range of applications make LIQUEEL® a truly versatile product.



Leg Ulcer on the Lower Leg

LOWER LEG **ULCUS CRURIS**

A typical wound on the lower leg, such as a leg ulcer, can be treated without having to remove individual foam segments.

The illustrative images show that LIQUEEL® can be applied very flexibly to flat, concave, and convex areas of the body.



Ulcer—covered with LIQUEEL®

L+F Medizinprodukte GmbH
Am Gewerbepark 1a
63594 Hasselroth
Germany [DE]

Phone: +49 6184 99 53 357
Fax: +49 6184 99 53 358
E-Mail: info@lf-med.de
Web: www.lf-medizinprodukte.de

