LF MED



CPK COMPRESS-2

2-LAYER COMPRESSION SYSTEM

26

2026./09

CPK COMPRESS-2


Two-layer compression system

CPK COMPRESS-2 is a two-layer compression system consisting of a soft, conformable padded bandage and an elastic, cohesive bandage, which together ensure uniform compression.

The padded bandage is soft, conforms to the skin, and is made of cotton, viscose, polyamide, and spandex. The cohesive bandage is skin-colored, made of 100% cotton, and coated on both sides with a latex-free cohesive adhesive. It is elastic, allowing the tension to be adjusted as needed and ensuring effective compression. Thanks to its cohesive and non-slip properties, the system stays in place for up to seven days.

How Does Compression Work

- Targeted external pressure constricts the dilated veins, allowing the venous valves to close again.
- The muscle-vein pump is activated, which increases the blood flow velocity and accelerates the return flow to the heart.
- Fluid buildup in the tissues (edema) is effectively reduced, and the healing process of chronic wounds is noticeably accelerated.

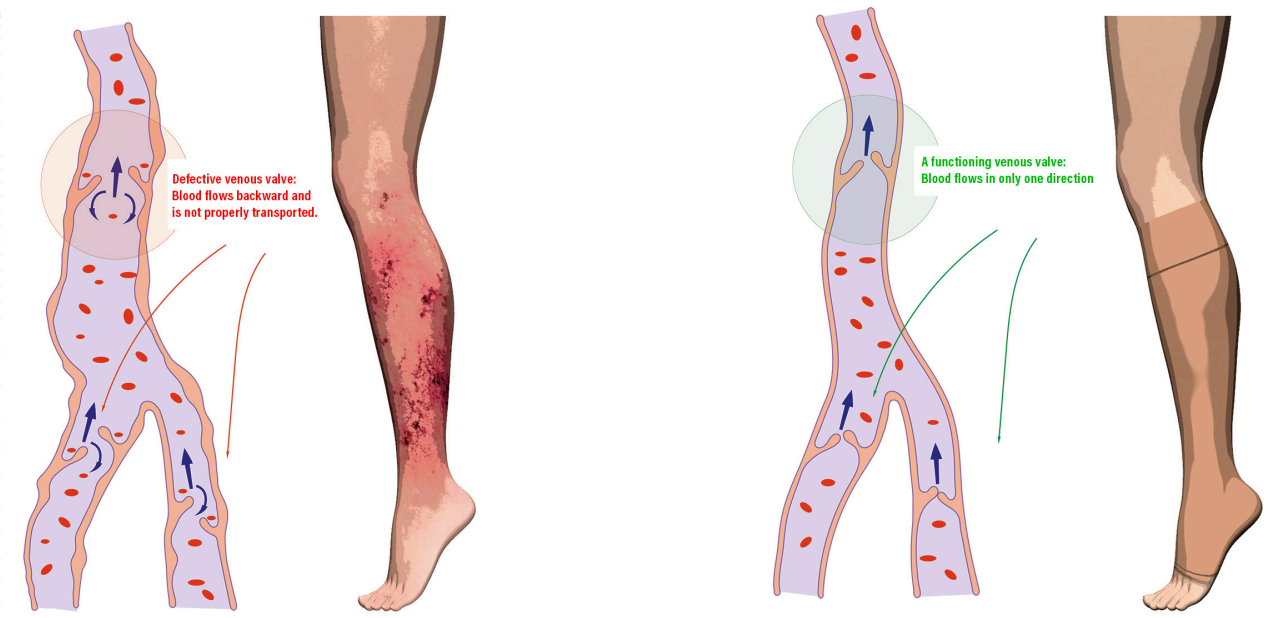

About 75% of all patients with venous and lymphatic vessel disorders can be successfully treated with compression therapy alone. Chronic leg ulcers (ulcus cruris) have been shown to heal significantly faster with consistent compression.

Performance

- ✓ *Precise pressure control of 20 or 40 mmHg*
- ✓ *Visual indicators for proper overlap*
- ✓ *High comfort thanks to latex-free padding*
- ✓ *Secure fit for up to 7 days*

CPK COMPRESS-2

Full-pressure and reduced-pressure compression systems

CPK COMPRESS-2 is used to treat venous leg ulcers, chronic venous edema, and other clinical conditions for which controlled compression is indicated. The system is available in various sizes and pressure levels, allowing its use to be individually tailored to patients' needs.

The reduced-pressure compression system CPK COMPRESS-2 Lite applies approximately 20 mmHg of pressure, while the full-pressure compression system CPK COMPRESS-2 reaches 40 mmHg.



Performance

- ✓ Precise pressure control of 20 or 40 mmHg
- ✓ Visual indicators for proper overlap
- ✓ High comfort thanks to latex-free padding
- ✓ Secure fit for up to 7 days


Indications

- Venous leg ulcers
- Lymphatic edema
- Other clinical conditions in which compression therapy is medically appropriate

	Size	Additional Information	Padding / Compression	Unit*	PZN**	REF
S	18 cm bis 25 cm	Full-Pressure System	10 cm x 4 m / 10 cm x 5.5 m	1 x 2	20362107	CPK-C-1
M	25 cm bis 32 cm	Full-Pressure System	10 cm x 5 m / 10 cm x 7.5 m	1 x 2	20362113	CPK-C-2
S	18 cm bis 25 cm	Pressure-Reduced System	10 cm x 4 m / 10 cm x 5.5 m	1 x 2	20368587	CPK-C-LT-1
M	25 cm bis 32 cm	Pressure-Reduced System	10 cm x 5 m / 10 cm x 7.5 m	1 x 2	20368601	CPK-C-LT-2

**Packaging unit; ** Pharmacy Central Number*

The Pharmazentralnummer (PZN) is a Germany-wide identification key for medicines and other pharmacy products.

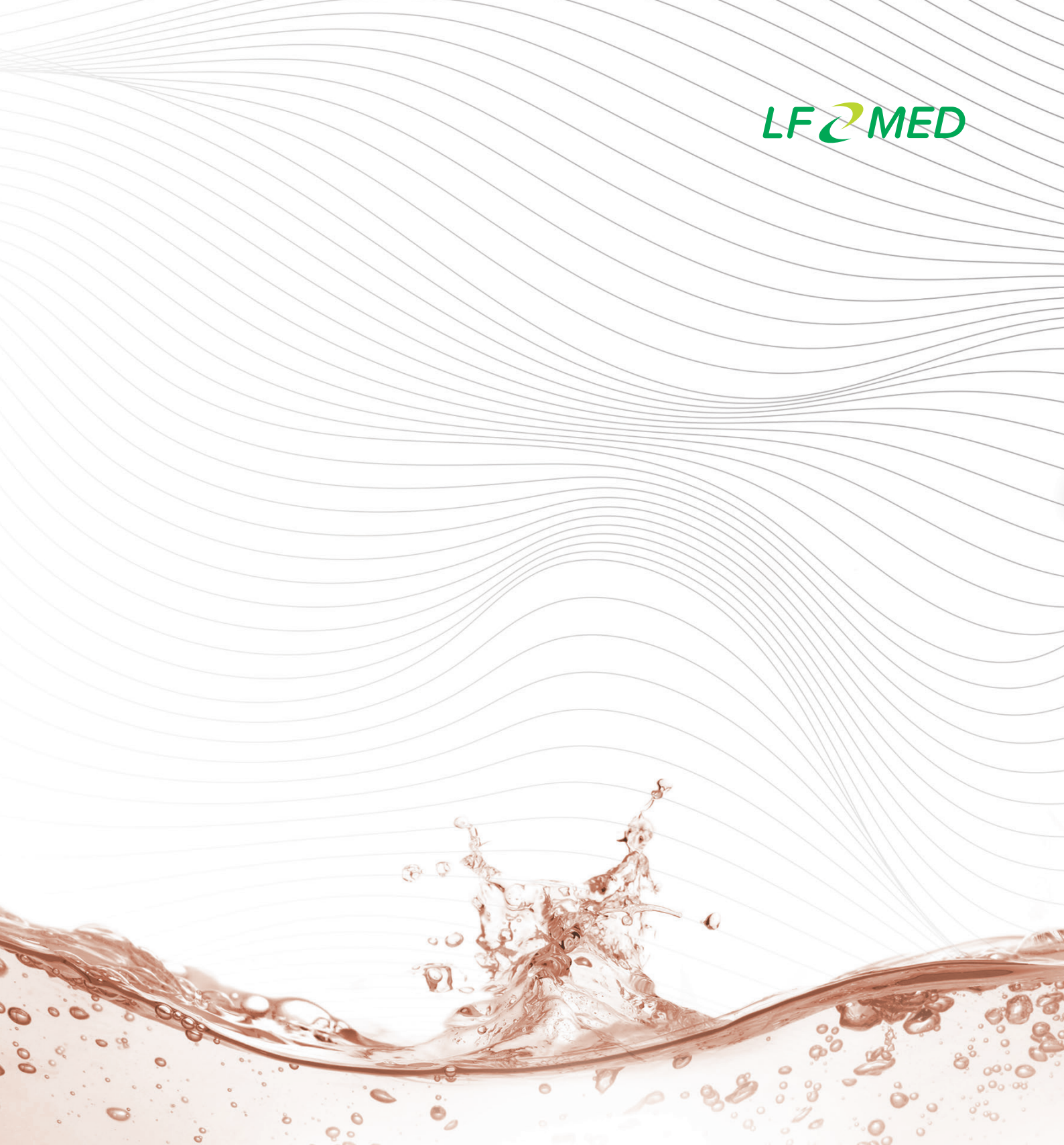

L+F Medizinprodukte GmbH
Am Gewerbepark 1a
63594 Hasselroth
Germany [DE]

Phone: +49 6184 99 53 357
Fax: +49 6184 99 53 358
E-Mail: info@lf-med.de
Web: www.lf-medizinprodukte.de